





Olink[®] Flex

Flexibility without compromise

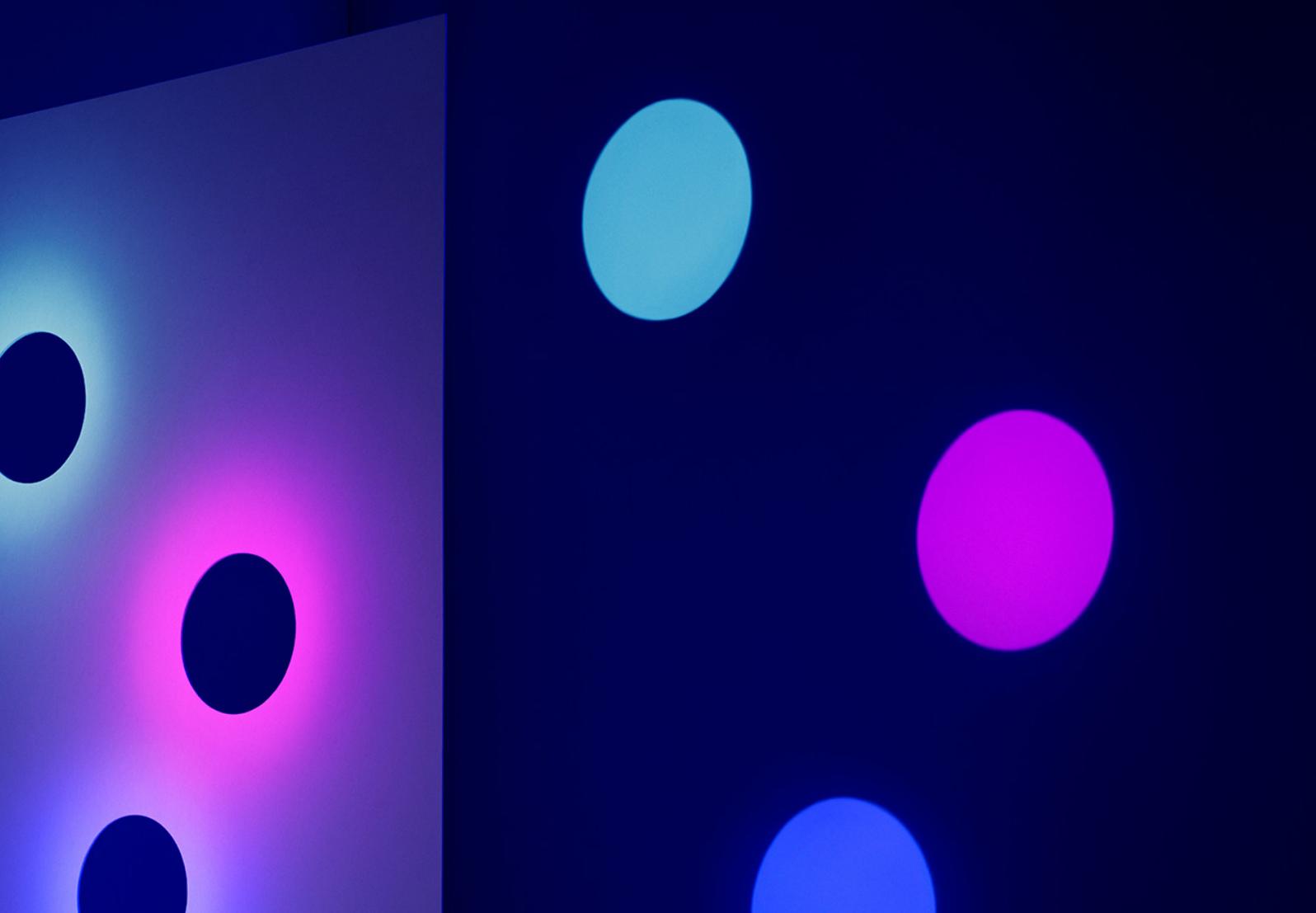

Olink® Flex — Inflammation research on your terms

Olink's highly flexible made-to-order product, Olink® Flex, enables researchers to select and combine targets for up to 21 human proteins in one biomarker panel. Pick and choose from over 200 inflammation-related human proteins in a broad mix-and-match library.

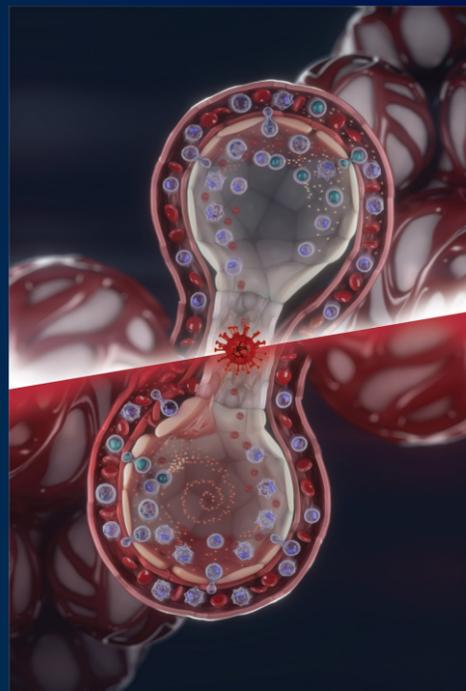

Inflammation and disease

Inflammation is activated in response to tissue damage, pathogens, toxins, and other factors. Persistent, chronic inflammation is associated with many diseases characterized by tissue and organ damage caused by the infiltration of primary inflammatory immune cells that release pro-inflammatory cytokines.

To help study inflammation across a broad spectrum of diseases, Olink offers the broadest and most comprehensive set of high-quality, multiplexed inflammation biomarker assays on the market. The list of available markers for Olink Flex can be found on the [Olink website](#).

How does it work?

Choose and combine 15-21 pre-validated assays from an extensive Olink library using our online panel design tool in Olink® Insight. Each assay maintains the quality, performance, specificity, and sensitivity that our other products are known for. Our internal controls remove the need to run replicate samples, and the data is reported in absolute (pg/mL) or relative (NPX) quantification.

Olink Flex gives you a streamlined workflow with robust performance. There is no need for time-consuming washing steps that could cause loss of material. Olink Flex runs on Olink® Signature Q100, a low maintenance, auto-calibrated instrument for data readout that facilitates unattended operation. Use the NPX Signature software for built-in QC and data analysis.



Uniquely scalable

The Olink platforms scale from our highplex product Olink® Explore 3072, to midplex Olink® Target 48 and Olink® Target 96 down to customized lowplex Olink® Focus and Olink® Flex.

All Olink platforms use our patented Proximity Extension Assay (PEA) technology, which provides exceptional specificity and sensitivity due to its dual recognition, DNA-coupled methodology. We offer the same high data quality regardless of plex-size.

High combinability

Select 15-21 protein biomarkers and analyze 40 samples simultaneously.

Convenient workflow

Less than 2 hours hands-on time and no washing steps.

Inflammation assay-focused

Biologically relevant inflammation-related proteins, selected in collaboration with experts.

Quantification method

Results reported in pg/mL (absolute quantification) or NPX (relative quantification).

Minimal sample volume

Use only 1 µL of blood to measure up to 21 proteins simultaneously.

Fast delivery

Receive your Olink Flex panel 3-5 weeks after placing the order.

What to expect when ordering Olink® Flex

Olink website

Visit the [Olink website](#) to find:

- Further product information
- Assay list
- Look at selected template panels
- Instruction video for the panel selection tool in Olink Insight
- Ordering process
- Info on show cases
- Validation data

Olink Insight

Visit the Flex panel builder and your own workspace at insight.olink.com to:

- Design and request your Olink Flex panel
- Select from template panels put together by experts
- Export and share panel designs with collaborators
- Request your selected panels

Olink Flex team

- Get contacted by an Olink expert for execution of your order
- Your panels are prepared and quality controlled

Delivery

- Kits, Panel Data File and a Certificate of Analysis (CoA) are delivered within 3-5 weeks
- The Panel Data file and CoA are delivered to your account in Olink Insight

www.olink.com

Olink products and services are For Research Use Only and not for Use in Diagnostic Procedures.
This document is not intended to convey any warranties, representations and/or recommendations of any kind, unless such warranties, representations and/or recommendations are explicitly stated.
Olink assumes no liability arising from a prospective reader's actions based on this document.
OLINK and the Olink logotype are trademarks registered, or pending registration, by Olink Proteomics AB.
© 2022 Olink Proteomics AB. All third-party trademarks are the property of their respective owners.
Olink products and assay methods are covered by several patents and patent applications <https://www.olink.com/patents/>.
Olink Proteomics, Dag Hammarskjölds väg 52B, SE-752 37 Uppsala, Sweden

1274, v1.0, 2022-11-09